



IAP4 Rec'd PCT TO 25 NOV 2005 PCT

IN THE UNITED STATES PATENT AND TRADEMARK OFFICE


Application No. : 10/537,027
Applicant : VAN DEN HEUVEL, ET AL.
Filed : MAY 31, 2005
Title : CLINICAL ASSISTANT FOR COCHLEAR IMPLANT CARE
Art Unit : TO BE ASSIGNED
Examiner : TO BE ASSIGNED
Atty Docket No. : COCH-0148-US1

Commissioner for Patents
P. O. Box 1450
Alexandria, VA 22313-1450
Sir:

The below-identified communication(s) is (are) submitted in the above-captioned application or proceeding:

- ☒ Certified copy of Australian Application No. 2002952675
- ☒ The Commissioner is hereby authorized to charge payment of any fees associated with this communication, including fees under 37 C.F.R. §§ 1.16 and 1.17 or credit any overpayment to **Deposit Account Number 10-0233-COCH-0148-US1**.

Respectfully submitted,

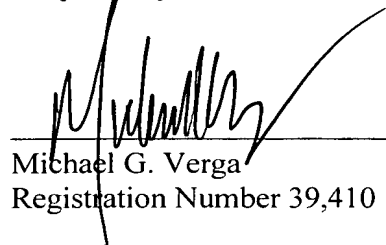

Michael G. Verga
Registration Number 39,410

JAGTIANI + GUTTAG
Democracy Square Business Center
10363-A Democracy Lane
Fairfax, Virginia 22030
(703) 591-2664

November 23, 2005



Australian Government

Patent Office
Canberra

I, JANENE PEISKER, MANAGER EXAMINATION SUPPORT AND SALES hereby certify that annexed is a true copy of the Provisional specification in connection with Application No. 2002952675 for a patent by COCHLEAR LIMITED as filed on 07 November 2002.

WITNESS my hand this
Twentieth day of October 2005

A handwritten signature in cursive script, appearing to read 'J. Peisker'.

JANENE PEISKER
MANAGER EXAMINATION SUPPORT
AND SALES



**CERTIFIED COPY OF
PRIORITY DOCUMENT**

AUSTRALIA
Patents Act 1990

Cochlear Limited

PROVISIONAL SPECIFICATION

Invention Title:

Clinical assistant for cochlear implant care

The invention is described in the following statement:

Technical Field

The present invention relates to a system for providing and delivering more efficient cochlear implant care.

5

Background of the Invention

Over recent times, hearing prosthesis, and in particular, cochlear implants have become more widespread in their use as the benefits become more widely realised by the hearing impaired.

10

Hearing loss can be due to many different causes. One type of hearing loss is conductive hearing loss which occurs when the normal mechanical pathways for sound to reach the hair cells in the cochlea are impeded, for example, by damage to the ossicles. Conductive hearing loss may often be helped by use of conventional hearing aids, which amplify sound so that acoustic information does reach the cochlea and the hair cells.

15

In many people who are profoundly deaf, however, the reason for deafness is sensorineural hearing loss. This type of hearing loss is due to the absence of, or destruction of, the hair cells in the cochlea which transduce acoustic signals into nerve impulses. These people are thus unable to derive suitable benefit from conventional hearing aid systems, because there is damage to or absence of the mechanism for nerve impulses to be generated from sound in the normal manner.

20

Cochlear implant systems have been developed for persons with sensorineural hearing loss which bypass the hair cells in the cochlea and directly deliver electrical stimulation to the auditory nerve fibres, thereby allowing the brain to perceive a hearing sensation resembling the natural hearing sensation normally delivered to the auditory nerve. US Patent 4532930, the contents of which are incorporated herein by reference, provides a description of one type of traditional cochlear implant system.

30

Typically, cochlear implant systems consist of a microphone for detecting acoustic signals and a sound processor for transforming the acoustic signals,

35

particularly speech, into patterns of electrical stimulation. The sound processor is typically worn externally, either behind the ear of the recipient or in a body worn pouch, and is programmed to meet the particular requirements of the recipient. The sound processor can include several different schemes for
5 processing the acoustic signal and transforming the signal into electrical stimuli, with these schemes well known in the art. The electrical stimuli is then transferred, together with a power signal, to an implanted receiver/stimulator unit positioned within the head of the user. Traditionally, this transcutaneous transmission occurs via the external transmitter coil which is positioned to
10 communicate with an implanted receiver coil provided with the receiver/stimulator unit. This communication serves two essential purposes; firstly to transcutaneously transmit the coded sound signal and secondly to provide power to the implanted receiver/stimulator unit. Conventionally, this link has been in the form of a radio frequency (RF) link, but other such links
15 have been proposed and implemented with varying degrees of success. The implanted receiver/stimulator unit traditionally includes a receiver coil that receives the coded signal and power from the external processor component, and a stimulator that processes the coded signal and outputs a stimulation signal to an intracochlea electrode assembly which applies the electrical
20 stimulation directly to the auditory nerve producing a hearing sensation corresponding to the original detected sound.

Due to the complex biophysical phenomena associated with the electrical excitation of neurons and the psychophysical phenomena regarding
25 the interpretation of neural activity by the auditory nervous system, not all recipients have been found to benefit from the same speech processing strategy. In this regard, it has been found that the quality and intelligibility of speech percepts evoked by a cochlear prosthesis may be improved in a given recipient by more specific manipulations of the electrical stimulus waveforms
30 tailored to that recipient.

In this regard, it is important that when a recipient is first implanted with a cochlear implant, in order for them to obtain the most benefit from the device, it needs to be adjusted/fitted to suit their specific needs. As the useful dynamic
35 range for electrical stimulation is relatively narrow and varies across recipients and electrodes, there is a need to individually tailor the characteristics of

electrical stimulation for each recipient. Simple psychophysical measurements establish the useful range for each electrode, and such parameters can be stored within the recipient's sound processor for continual use. This procedure is often referred to as "mapping" and is the term given to the process of measuring and controlling the amount of electrical current delivered to the cochlea. It is this process that ensures stimulation from the implant provides a recipient with comfortable and useful auditory percepts, and is essential in ensuring that the recipient receives maximum benefit from the cochlear implant.

As the implant system is designed to present acoustic information, in particular speech, to a recipient in a usable form, the initial aim of the mapping process is to optimise the information provided for a particular recipient. The "mapping" process is a key part of the post-operative management of all cochlear implant recipients and occupies a significant proportion of the post-operative clinical time.

As well as ensuring that the implant is fitted to the recipient immediately following implantation, it is also important to provide continual, ongoing management of the device. Typically, the initial fitting session is aimed at providing the recipient with the ability to experience sounds produced by the implant, so that they can become accustomed to the type of sounds experienced. After this stage, it is important that the sound information is optimised to ensure that the benefit obtained from the cochlear implant is maximised. This is particularly important during the first three months following initial fitting as well as over the lifetime of the recipient at regular six month intervals or when the recipient requests an update or reports an issue.

Lifelong after-care for the recipient is essential in ensuring the continuing success of the implant and that the recipient obtains maximum benefit from the device. Such after-care includes; continual or regular monitoring of the operation of the device and comparison of current performance against historical records for that device; continual checking of the stability of the technical and physiological operation of the device; updating the device as new technology becomes available; ongoing technical service of the external parts of the system; distribution of spare parts and replacement batteries; as well as ongoing recipient counselling. Such follow-up sessions typically require

standard checks of the recipient's threshold and comfort levels and the integrity of the overall system including speech processor and microphone.

During after-care, it is important that any problems are quickly and efficiently detected and remedied. If there are areas where the recipient is experiencing difficulty, or where tests suggest a problem, such as speech perception, it is important that this is identified and investigated as the growth and development of the recipient's communication skills may become affected. This is particularly important in relation to infants implanted with cochlear implant systems to ensure that they can develop speech and communication skills in their early development phase.

Studies have shown that in some clinics that provide after-care to cochlear implant recipients, this function alone can take up to 30% of the clinic's total cochlear implant program workload. As the after-care function is provided mostly by audiologists, this maximisation of the audiologist's time directly competes with the time available for the audiologist to identify and work with potential new recipients who may benefit from a cochlear implant.

Typically, much of the after-care requires dealing with relatively minor issues, such as hardware defects and/or use of accessories or replacements. Such after-care does not necessarily require the attention of an audiologist as their time is best spent on adjusting a recipient's map or speech processing settings, in response to a change in hearing status.

Further, many recipients live reasonable distances from support clinics, which are mainly located in larger cities, and as such they may not visit a clinic as often as a recipient living relatively close to a clinic. Such recipients may only really contact the clinic if they have a specific problem which greatly affects the performance of their implant, rather than routinely. This may result in such recipients not obtaining maximum benefit from their implant.

The present invention therefore aims to address these problems.

Any discussion of documents, acts, materials, devices, articles or the like which has been included in the present specification is solely for the purpose of

providing a context for the present invention. It is not to be taken as an admission that any or all of these matters form part of the prior art base or were common general knowledge in the field relevant to the present invention as it existed before the priority date of each claim of this application.

5

Summary of the Invention

Throughout this specification the word "comprise", or variations such as "comprises" or "comprising", will be understood to imply the inclusion of a
10 stated element, integer or step, or group of elements, integers or steps, but not the exclusion of any other element, integer or step, or group of elements, integers or steps.

The present invention preferably provides a system that allows for
15 quicker and more efficient follow-up care of recipients of cochlear implants that enables a full, standard and consistent service to be provided that increases the productivity of the audiologists during the clinical time; increases the quality of the follow-up service; standardises the results of follow-up tests; and reduces the number of erroneously diagnosed problems.

20

The present invention also preferably provides to the recipients of cochlear implants increased quality of follow-up tests; a sense of control and participation in the monitoring of the performance of their device; and a system which allows for follow-up services to be done in a more friendlier environment,
25 such as in the recipient's home or in a hearing centre.

It is considered that by providing an after-care system that is more efficient and user friendly as well as less audiologist-dependent, a more cost effective and attractive cochlear implant management process is able to be
30 implemented.

The present invention also preferably provides a system for use by recipients in follow-up visits that provides a relatively easy to use method of conducting such tests.

35

In a first aspect, the present invention is a system for performing one or more tests on a hearing prosthesis that is usable at least in part by the recipient of the prosthesis, the system comprising:

5 a computer that processes software instructions and outputs signals in response to said instructions;

a prosthesis interface means that provides transfer of signals from said computer to a prosthesis and/or from the prosthesis to the computer; and

10 a recipient computer interface device that allows the recipient of the prosthesis to at least partially control at least some aspect of the tests performed on a prosthesis that is interfaced with the computer.

In this aspect, the computer can in one embodiment be adapted to stand alone and perform said one or more tests on the hearing prosthesis. In this case, the computer may be housed in a clinic of the clinician and then used by
15 the recipient when visiting their clinician. The results of the tests can be assessed immediately by the clinician or at some later time that is perhaps more convenient for the clinician.

In another embodiment, the computer can be housed at a location
20 remote from the clinician. In a preferred embodiment, the computer is housed in the recipient's home or at a location that is relatively more convenient for the recipient to access than the clinic. In this case, the results of the tests can be stored in a storage means in the hearing prosthesis or can be stored in a suitable storage means used in conjunction with the computer. For example,
25 the results can be saved on a portable storage means, such as a disc, and then provided to a clinician responsible for managing the after-care of the recipient of the prosthesis. The stored results may be retained by the recipient and only provided to the clinician on the recipient's next visit to the clinician or can be posted to the clinician using normal mail services.

30

In a further embodiment, the clinician can set up the configuration for one or more tests in the hearing prosthesis. When the recipient decides or is instructed to test their prosthesis, the prosthesis can be interfaced with the computer with the computer being adapted to read the test configuration from
35 the hearing prosthesis.

In another embodiment, the computer is adapted to deliver the results electronically to a computer used by the clinician. Such electronic transmission is preferably provided by the Internet but other suitable modes of transmission can be envisaged. In this embodiment, the clinician's computer is preferably
5 remote from the first computer that is interfaced to the hearing prosthesis. The clinician's computer can preferably store the results delivered from the first computer for review and assessment by the clinician at a convenient time.

The clinician's computer can further be adapted to interact with the first
10 computer. In this embodiment, the clinician's computer can provide a set of instructions to the first computer instructing the first computer to perform one or more tests on the hearing prosthesis. The instructions of the clinician's computer can be based on or modified depending on the results received from the first computer. In one embodiment, the first computer can be essentially a
15 dumb terminal and act mainly as a means of delivering the instructions from the clinician's computer. In a more preferred embodiment, the first computer simply requires a trigger signal from the clinician's computer which then activates appropriate software already loaded on the first computer which results in the first computer performing one or more desired tests on the
20 prosthesis.

In a second aspect, the present invention is a system for performing one or more tests on a hearing prosthesis comprising:

a first computer that processes a set of software instructions and outputs
25 signals in response to said instructions;

a prosthesis interface means that provides transfer of signals from said first computer to a prosthesis and/or from the prosthesis to the computer; and

a second computer that provides said set of instructions to said first computer to control at least some aspects of the tests performed on a
30 prosthesis that is interfaced with said first computer.

In this aspect, the first computer can be housed at a location remote from the clinician. In a preferred embodiment, the first computer is housed in the recipient's home or at a location that is relatively more convenient for the
35 recipient to access than the clinic. The first computer is preferably adapted to deliver the results electronically to the second computer that is preferably used

by the clinician. Such electronic transmission is preferably provided by the Internet but other suitable modes of transmission can be envisaged. The second computer can preferably store the results delivered from the first computer for review and assessment by the clinician at a convenient time.

5

The second computer can further be adapted to interact with the first computer. In this embodiment, the second computer can provide a set of instructions to the first computer instructing the first computer to perform one or a series of tests on the hearing prosthesis. The instructions of the second
10 computer can be based on or modified depending on the results received from the first computer. In one embodiment, the first computer can be essentially a dumb terminal and act mainly as a means of delivering the instructions received from the second computer. In a more preferred embodiment, the first computer simply requires a trigger signal from the second computer which then
15 activates appropriate software already loaded on the first computer which results in the first computer performing one or more desired tests on the prosthesis.

In the second aspect, the system further preferably comprises a
20 recipient-first computer interface device that allows the recipient of the prosthesis to at least partially control at least some aspect of the tests performed on a prosthesis that is interfaced with the computer.

In both aspects, the recipient interface device preferably comprises a
25 monitor, a keyboard, a keypad, and/or a pointing device, such as a mouse or stylus. For example, the monitor preferably displays a number of messages to the recipient. These messages can include instructions as to how to perform the test. The monitor also preferably provides messages advising the recipient how to commence and/or stop a test. A test can typically be started or stopped
30 by pressing any key of the computer's keyboard or by clicking on an icon displayed on the monitor using the pointing device.

In either aspect, the hearing prosthesis is preferably a cochlear implant. For the purposes of the description herein, the systems will be described more
35 fully in relation to their preferred use in testing cochlear implants. It will,

however, be appreciated that the systems could be used to perform tests on hearing aids.

The second computer is preferably loaded with appropriate software that
5 enables the clinician operating this computer to configure an appropriate set of instructions that can be sent to the recipient's computer over the appropriate network connection. To allow this, the second computer is also preferably provided with an interface that enables the clinician to configure a series of tests to measure desired aspects of the recipient's cochlear implant system.
10 The interface of the second computer again preferably comprises a display provided on a monitor that can be manipulated or controlled by the clinician using a keypad, keyboard and/or a pointing device, such as a mouse.

The prosthesis interface means can comprise a socket capable of being
15 connected via a cable to the external speech processor of the cochlear implant of the recipient. It is envisaged that the speech processor would during most if not all tests be mounted so as to be in signal communication with the implanted stimulator unit of the implant, thereby providing complete access of the system to the cochlear implant.

20

In a further embodiment, the first computer can also be provided with a custom input device for use by the recipient. Such a device may include a series of buttons and/or lights that require manipulation by the recipient during testing of the cochlear implant.

25

Upon delivery of the set of instructions to the first computer, the recipient can preferably at least partially control the instructions by reading the messages displayed on the monitor and responding appropriately. The monitor, for example, may display a message asking the recipient to perform a
30 series of tasks. The computer then collects relevant information resulting from the completion of such tasks.

Once the one or more tests have been completed, a set of results collected during the execution of the tests is then preferably sent back to the
35 clinician's computer for evaluation by the clinician. In this embodiment, the clinician's computer can perform an appropriate evaluation of the data received

and assess it against predefined criteria set by the clinician. The results of the assessment can then be conveyed to the clinician who is then in a position to make an informed decision on what, if any, subsequent action can be made.

5 The step of evaluating the received results can have a number of outcomes, depending upon the results of the evaluation. Should the evaluation show that the performed tests were passed, no further action will be required by the recipient and the recipient would normally be informed of this. However, should the results of the tests indicate a problem with the cochlear implant that
10 requires corrective action, the recipient would be informed of the need to visit the clinic/hospital and an appointment could be made and confirmed by the recipient for such a session. With the present invention, it is also possible that should the test results be inconclusive, the clinician could configure further tests to ascertain the function of the cochlear implant device.

15

 The monitor of the clinician's computer can output a first display that assists the clinician in configuring the appropriate set of instructions for delivery to the recipient's computer. This first display allows the clinician to enter the details of a recipient or select a recipient from a database of recipients and also
20 to view the recipient's detailed data and history of past tests performed, to establish whether there are any specific issues that require testing.

 Once the clinician has selected the recipient and considered their history, they must then configure appropriate tests to be undertaken to assess the
25 performance and the operation of the cochlear implant. In a further embodiment, the displays provided on the clinician's computer and/or the recipient's computer can display messages in the preferred language suitable for the specific recipient. The software preferably allows ready changing of the language of the messages displayed by one or both computers.

30

 Once the clinician has accessed the first display, the clinician preferably has three options. They can create a new test configuration, open a template test configuration, or open a recipient test configuration. This decision will be largely dependant upon the knowledge the clinician possesses of the
35 recipient's cochlear implant system.

It is considered that in most instances, the clinician will open a template test configuration which preferably provides a standard template for configuring the tests to suit the chosen recipient.

- 5 Should the clinician have a previous test template stored that has been configured for the particular recipient, the recipient test configuration would normally be chosen. Alternatively, should the clinician wish to construct a new test configuration different to that of the template, the option to create a new test configuration will be selected.

10

Following this step, the tests need to be configured and this can be done by manipulating the template displayed on the clinician's computer. By doing this, the clinician is able to run one or more standard tests checking a number of features of the cochlear implant system. These could be:

15

(i) a standard system integrity test to ensure that all the components of the system are operating correctly;

(ii) a neural response threshold test to compare current measured
20 neural response thresholds with previously measured thresholds;

(iii) threshold (T) stimulation level and comfortable (C) stimulation level test to determine whether the correct dynamic range of stimulation is set for each electrode; and

25

(iv) tests to check loudness balancing, pitch ranking and phoneme discrimination to give the clinician an indication of the effectiveness of the implant.

30 The clinician can also select which stored maps are to be used in the testing.

Alternatively, in this step the clinician is able to select an advanced configuration to customise the standard tests to suit the specific needs of the
35 recipient. In this regard, the clinician could modify the input device used, the specific electrodes used in the threshold level tests, the acceptable difference

allowable between measured results and stored data, the order of tests to be performed, and the samples taken for the speech tests. Once the clinician has selected such a configuration, the clinician can then alter specific aspects of the tests to suit the needs of the recipient.

5

Once all the tests have been configured and set for the specific recipient, the clinician can save the configuration for remote execution, to be sent to the recipient's computer.

10 Should the recipient be present with the clinician at the time of configuration, the clinician can also start testing the recipient by starting the recipient's interface upon completion of the configuration session. In this instance, the clinician can check whether the speech processor (SP) is connected to the computer by running a speech processor connection test, to
15 ensure that the test collects desired results.

Following delivery of the configured tests at the recipient's computer, the recipient runs the test via an appropriate interface, such as has been defined herein.

20

The recipient's computer takes the recipient through a series of tests in sequence, as determined by the clinician, collecting the data and results of each test. After each test is performed, the measured results are merged into a results file where this file is then used by the clinician's computer following the
25 completion of the recipient session. This file can be exported to the clinician in a number of ways, for example the file can be sent via the Internet or LAN connection to the recipient where it is received by the clinician application and acted upon. Or it can be stored in the recipient's speech processor and accessed from the speech processor at the recipient's next visit to the clinician.

30

An example configuration file is shown below:

```
<?xml version="1.0" encoding="utf-8" ?>
```

```
<Configuration>
```

```
35   <PatientID>AEB182FC-F435-4f7e-85AE-D0D262464D86</PatientID>
      <TestID>B7316444-052D-464c-BC68-DC6474B822AA</TestID>
      <TestItems>
```

```

5      <Test1>
        <Name>Impedance Test</Name>
        <Channels>
          <Channel>1</Channel>
          <Channel>5</Channel>
        </Channels>
        <Modes>
          <Mode>CG</Mode>
          <Mode>MP1</Mode>
        </Modes>
10     </Test1>
      <Test2>
        <Name>Threshold Level Measurement</Name>
        <Map>1</Map>
15     <Channels>
          <Channel>1</Channel>
          <Channel>5</Channel>
          <Channel>10</Channel>
        </Channels>
20     </Test2>
    </TestItems>
  </Configuration>

```

This example configures an impedance test (test 1) and a threshold level
 25 measurement test (test 2) on the recipient application. As can be seen, the
 impedance test is carried out on channels 1 and 5 of the recipient's cochlear
 implant in both common ground (CG) and MP1 mode. The threshold level
 measurement test is carried out on channels 1, 5 and 10. Each of these
 parameters are set by the clinician when configuring the test.

30

Based upon this configuration file shown above, an example results file
 generated is also shown below.

```

    <?xml version="1.0" encoding="utf-8" ?>
35  <Result>
    <PatientID>AEB182FC-F435-4f7e-85AE-D0D262464D86</PatientID>
    <TestID>B7316444-052D-464c-BC68-DC6474B822AA</TestID>
    <TestItems>
      <Test1>
40      <Name>Impedance Test</Name>
        <Impedance>
          <Mode>CG</Mode>
          <Channel>1</Channel>
          <Value>10500</Value>

```



```

5      </Impedance>
      <Impedance>
        <Mode>CG</Mode>
        <Channel>5</Channel>
        <Value>10560</Value>
      </Impedance>
      <Impedance>
        <Mode>MP1</Mode>
        <Channel>1</Channel>
        <Value>10900</Value>
      </Impedance>
      <Impedance>
        <Mode>MP1</Mode>
        <Channel>5</Channel>
        <Value>10800</Value>
      </Impedance>
      </Test1>
      <Test2>
        <Name>Threshold Level Measurement</Name>
        <Level>
          <Channel>1</Channel>
          <Value>105</Value>
        </Level>
        <Level>
          <Channel>5</Channel>
          <Value>115</Value>
        </Level>
        <Level>
          <Channel>10</Channel>
          <Value>106</Value>
        </Level>
      </Test2>
    </TestItems>
  </Result>
35

```

As is shown above, the results of the configured tests, test 1 and test 2, are merged into this results file.

Following collection of the data for each test, the results data is preferably then sent back to the clinician for evaluation and decision on what course of action is required.

Once the recipient has received the configured test from the clinician, the recipient's computer preferably displays a personal welcome page to the

recipient wherein the recipient can be satisfied that the test is specifically designed for their own needs. The recipient can control commencement of the test by performing an action requested by the computer, for example by clicking the word "Next" using a mouse pointer. On doing this, the recipient is

5 preferably informed that the system will perform a system integrity check, to check that each aspect of the system is working correctly. The recipient is warned before any testing and in order for the test to be performed that recipient must select to proceed with the test by performing an action requested by the computer, such as clicking on the word "Go" using a mouse pointer.

10

Should the recipient experience any difficulties or discomfort during testing, they have the ability to stop the testing at any time by pressing any button or "Stop". This feature is similar to what would occur in a regular clinical visit, where the recipient could indicate to the clinician any issues and the

15 clinician would then stop the testing.

Once the recipient has started the system integrity test, the recipient is informed by a message on the monitor that the test is being performed and is provided with an indication of the progress of the test. Upon completion of the

20 test, the recipient is informed by a message on the monitor that the test has been completed and is asked to proceed to the next stage in the process. Again, the recipient preferably has control of commencement of the next stage of the test.

25 At the completion of the system integrity test, the recipient application preferably undergoes a basic test that checks whether the system integrity test has found any severe system integrity problems. Should such problems be found, the application would inform the recipient of this and suggest that the recipient contact the clinician for further attention. This would then normally

30 finish the testing process, as the presence of system integrity problems would typically reduce the benefit of obtaining any further test results.

Following the successful completion of the system integrity test, the recipient is then preferably informed of the next test in the process, such as a

35 psychophysics test. Upon commencement of this test, the recipient is given instructions on how the test will be conducted and direction as to how input is

required from the recipient. Typically, the input will be in the form of using a mouse to activate the appropriate icon on the computer screen. For example, the recipient can be asked to press an icon "YES" when they hear a tone and press an icon "NO" when they do not hear anything. Once this test is
5 complete, the recipient is preferably informed of this and asked to proceed to the next step in the process.

When the recipient has completed all tests configured for them, they are preferably informed of this and the application is closed. The recipient
10 preferably clicks an icon "End" to close the programme on their computer.

As discussed, upon completion of the configured tests, the results are preferably sent to the clinician for evaluation. This is done by the clinician application which accesses the results file and tests the results in relation to
15 previous or desired results.

Upon an evaluation of the results as discussed above for all the tests undertaken, the clinician is preferably provided with an interface on the clinician's computer providing an overview of the results. The clinician can then
20 decide to view a detailed report of the test. Such a detailed report can show the measured values obtained in the test against the mapped values, or previously measured values in both a data form and a graphical form.

Based upon the results, the clinician can decide to proceed in a number
25 of ways. Obviously, if the process indicated that all tests were passed, no specific action would be required and the clinician would inform the recipient of this and ask the recipient to repeat the test in, for example, 6 months time. Should there be some inconsistencies in the results, the clinician may wish to see the recipient personally or send the recipient further tests to undertake to
30 determine whether the inconsistency is merely an error in recordings or is a concern that requires further monitoring.

According to a still further aspect, the present invention is a method of testing a hearing prosthesis using a system as defined herein.

As can be appreciated, the present invention provides a system which enables faster, more user friendly and less skill-intensive follow-up care and management of a cochlear implant. This system allows for personal expert care from clinicians and audiologists without the need for onerous clinician visits, and allows the recipient to interact with the clinician from their own home or school.

Brief Description of the Drawings

By way of example only, preferred embodiments of the invention are now described with reference to the accompanying drawings, in which:

Fig. 1 is a view of one embodiment of the present system;

Fig. 2 is a flow diagram representing the overall process of the present invention;

Fig. 3 is a flow diagram representing the process associated with the clinical application of the present invention;

Fig. 4 is a representation of an embodiment of a clinical interface of the present invention;

Fig. 5 is an embodiment of an advanced clinical interface of the present invention;

Fig. 6 is a flow diagram of one embodiment of the recipient application of the present invention;

Figs. 7-14 depict an example of one embodiment of a test as performed by the recipient;

Fig. 15 shows an embodiment of an evaluation of a test process according to the present invention; and

Fig. 16 shows an embodiment of an advanced evaluation report generated as a result of the present invention.

Best Mode for Carrying Out the Invention

5

The following description is a preferred embodiment of the present invention and is not to be taken in a limiting sense but is shown purely for the purpose of describing the general principles of the invention.

10 The invention is generally embodied on a computer system including one or more computers, such as the set-up as shown in Fig. 1. In the depicted embodiment, two stand alone personal computers are provided, with one computer 2 being for use by a clinician and the other computer 4 being for use by a recipient. These computers can be located in positions remote from each
15 other and can be connected via an appropriate network connection, such as a LAN (Local Area Network) or via the Internet. In this regard, the computers can be located in the same room, such as a clinician's room; in the same building, such as in the same clinic or hospital; or in separate buildings or locations, with, for example, computer 2 being in a clinic and computer 4 being at the
20 recipient's home, school or workplace.

While two computers are depicted, it will be appreciated that it is also possible for the present invention to be performed on one computer that is used by both the clinician and the recipient at different times. In another
25 embodiment, the present invention could be performed using a single computer. In a further embodiment, the clinician can set up the configuration for one or more tests in the hearing prosthesis. When the recipient decides or is instructed to test their prosthesis, the prosthesis can be interfaced with a computer that is adapted to read the test configuration from the hearing
30 prosthesis. This computer can be stand-alone computer and can belong to the recipient, the clinician or a third party.

As is shown in Fig. 1, the clinician's computer 2 is a general stand-alone personal computer having a monitor, keyboard and mouse. This computer is
35 loaded with appropriate software according to the present invention to enable the clinician operating the computer to configure an appropriate recipient

application 6 that can be sent to the recipient's computer 4 over the appropriate network connection. This is done by providing the clinician with an interface that enables them to configure a series of tests to measure desired aspects of the recipient's cochlear implant system. This will be described in more detail below. Once the clinician has completed configuring the desired tests the application 6 is then sent to the recipient's computer 4.

The recipient's computer, like the clinician's computer 2, is a general stand-alone personal computer having a monitor, keyboard and mouse as well as connections to communicate with the recipient's cochlear implant system. This could be in the form of a simple socket capable of being connected to the recipient's external speech processor, which is in turn connected to the implanted stimulator unit, thereby providing complete access of the system to the cochlear implant.

As depicted, the computer 4 can also be provided with a custom input device 8 which may include a series of buttons and/or lights that require manipulation by the recipient during testing of the cochlear implant.

Upon delivery of a recipient application 6, the recipient can then run the application on the computer 4 and follow the instructions of the application. In this regard, the computer 4 is provided with application software including a recipient interface that receives the recipient application 6 and performs the required tasks and communicates this to the recipient during execution. This interface may ask the recipient to perform a series of tasks and collects the relevant information resulting from the completion of such tasks.

Once the recipient application 6 has been completed, the result data 9 collected during the execution of the recipient application 6 is then sent back to the clinician's computer 2 for evaluation by the clinician. In this regard, the clinician's computer 2 performs an appropriate evaluation of the data received and assesses it against predefined criteria set by the clinician. The results of the assessment is then conveyed to the clinician who is then in a position to make an informed decision on what, if any, subsequent action can be made.

This process is shown in Fig. 2. As is shown, the step of evaluating the received results can have a number of outcomes, depending upon the results of the evaluation. Should the evaluation show that the performed tests were passed, no further action will be required by the recipient and the recipient would be informed of this. However, should the results of the tests indicate a problem with the cochlear implant that requires corrective action the recipient would be informed of the need to visit the clinic/hospital and an appointment could be made and confirmed by the recipient for such a session. With the present invention, it is also possible that should the test results be inconclusive, the clinician could configure further tests to ascertain the function of the cochlear implant device.

The step of configuring the recipient application is now described in more detail with reference to Fig. 3. In this regard, an interface such as that shown in Fig. 4 is displayed on the monitor of the clinician's computer 2 to assist the clinician in configuring the appropriate application for the recipient's cochlear implant. This interface allows the clinician to select a recipient from a database of recipients and also to view the recipient's detailed data and history of past tests performed, to establish whether there are any specific issues that require testing.

Once the clinician has selected the recipient and considered their history, they must then configure appropriate tests to be undertaken to assess the performance and the operation of the cochlear implant. As shown in Fig. 4, the present invention can be used universally with the clinician being able to select the preferred language suitable for the specific recipient.

Once the clinician has accessed the display depicted in Fig. 4, the clinician has three options. They can create a new test configuration, open a template test configuration, or open a recipient test configuration. This decision is largely dependant upon the knowledge the clinician possesses of the recipient's cochlear implant system.

It is considered that in most instances, the clinician will open a template test configuration such as that shown in Fig. 4, which provides a standard template for configuring the tests to suit the chosen recipient. Should the

clinician have a previous test template stored that has been configured for the particular recipient, the recipient test configuration would normally be chosen. Alternatively, should the clinician wish to construct a new test configuration different to that of the template, the option to create a new test configuration will
 5 be selected.

Following this step, the tests need to be configured and this can be done by manipulating the template as shown in Fig. 4. By doing this, the clinician is able to run standard tests checking a number of features of the cochlear
 10 implant system. These could be:

- (i) a standard system integrity test to ensure that all the components of the system are operating correctly;
- (ii) a neural response threshold test to compare current measured
 15 neural response thresholds with previously measured thresholds;
- (iii) threshold stimulation level and comfortable stimulation level test to determine whether the correct dynamic range of stimulation is set for each electrode; and
- (iv) tests to check loudness balancing, pitch ranking and phoneme
 20 discrimination to give the clinician an indication of the effectiveness of the implant.

As is shown, the clinician can also select which stored maps are to be used in the testing.

25

Alternatively, in this step the clinician is able to select an advanced configuration to customise the standard tests to suit the specific needs of the recipient. In this regard, the clinician could modify the input device used, the specific electrodes used in the threshold level tests, the acceptable difference
 30 allowable between measured results and stored data, the order of tests to be performed, and the samples taken for the speech tests. Once the clinician has selected the advanced configuration, as shown in Fig. 5, the clinician can then alter specific aspects of the tests to suit the needs of the recipient.

Once all the tests have been configured and set for the specific recipient, the clinician can save the configuration for remote execution, to be sent to the recipient's computer 4.

- 5 As is also shown in Fig. 4, should the recipient be present with the clinician at the time of configuration, the clinician can also start testing the recipient by starting the recipient's interface upon completion of the configuration session. In this instance, the clinician can check whether the speech processor (SP) is connected to the computer 4 by running the SP
10 connection test, to ensure that the test collects desired results.

Figure 6 depicts the process associated with the recipient receiving and performing the configured test discussed in Figs. 4 and 5. As is shown, following delivery of the configured tests at the computer 4, the recipient runs
15 the test via an appropriate interface.

As is shown, the interface then takes the user through a series of tests in sequence, as determined by the clinician, collecting the data and results of each test. Tests 1-5 as shown in Fig. 6 each configure a test module in the
20 recipient application. After each test is performed, the measured results are merged into a results file where this file is then used by the clinician application following the completion of the recipient session. This file can be exported to the clinician in a number of ways, for example the file can be sent via the Internet or LAN connection to the recipient where it is received by the clinician
25 application and acted upon. Or it can be stored in the recipient's speech processor and accessed from the speech processor at the recipient's next visit to the clinician.

An example configuration file is shown below:

30 <?xml version="1.0" encoding="utf-8" ?>
<Configuration>
 <PatientID>AEB182FC-F435-4f7e-85AE-D0D262464D86</PatientID>
 <TestID>B7316444-052D-464c-BC68-DC6474B822AA</TestID>
35 <TestItems>
 <Test1>
 <Name>Impedance Test</Name>
 <Channels>

```

        <Channel>1</Channel>
        <Channel>5</Channel>
    </Channels>
    <Modes>
5        <Mode>CG</Mode>
        <Mode>MP1</Mode>
    </Modes>
</Test1>
<Test2>
10    <Name>Threshold Level Measurement</Name>
    <Map>1</Map>
    <Channels>
        <Channel>1</Channel>
        <Channel>5</Channel>
15        <Channel>10</Channel>
    </Channels>
</Test2>
</TestItems>
</Configuration>
20

```

This example configures an impedance test (test 1) and a threshold level measurement test (test 2) on the recipient application. As can be seen, the impedance test is carried out on channels 1 and 5 of the recipient's cochlear implant in both common ground (CG) and MP1 mode. The threshold level measurement test is carried out on channels 1, 5 and 10. Each of these parameters are set by the clinician when configuring the test as discussed in relation to Figs. 3, 4 and 5.

Based upon this configuration file shown above, an example results file generated is also shown below.

```

<?xml version="1.0" encoding="utf-8" ?>
<Result>
35    <PatientID>AEB182FC-F435-4f7e-85AE-D0D262464D86</PatientID>
    <TestID>B7316444-052D-464c-BC68-DC6474B822AA</TestID>
    <TestItems>
        <Test1>
            <Name>Impedance Test</Name>
            <Impedance>
40                <Mode>CG</Mode>
                <Channel>1</Channel>
                <Value>10500</Value>
            </Impedance>

```

```

5      <Impedance>
      <Mode>CG</Mode>
      <Channel>5</Channel>
      <Value>10560</Value>
      </Impedance>
      <Impedance>
      <Mode>MP1</Mode>
      <Channel>1</Channel>
      <Value>10900</Value>
10     </Impedance>
      <Impedance>
      <Mode>MP1</Mode>
      <Channel>5</Channel>
      <Value>10800</Value>
15     </Impedance>
    </Test1>
    <Test2>
      <Name>Threshold Level Measurement</Name>
      <Level>
      <Channel>1</Channel>
      <Value>105</Value>
20     </Level>
      <Level>
      <Channel>5</Channel>
      <Value>115</Value>
25     </Level>
      <Level>
      <Channel>10</Channel>
      <Value>106</Value>
30     </Level>
    </Test2>
  </TestItems>
</Result>

```

35 As is shown above, the results of the configured tests, test 1 and test 2, are merged into this results file.

Following collection of the data for each test, the results data is then sent back to the clinician for evaluation and decision on what course of action is
40 required.

An example recipient interface that would appear on the recipient's computer 4 running the test configuration prepared in Fig. 4 will now be described in relation to Figs. 7 - 14.

As shown in Fig. 7, once the recipient has received the configured test from the clinician, the interface runs a personal welcome page to the recipient wherein the recipient can be satisfied that the test is specifically designed for their own needs. Once the recipient proceeds by clicking "Next", the recipient is informed, as depicted in Fig. 8, that the system will perform a system integrity check, to check that each aspect of the system is working correctly. The recipient is warned before any testing and in order for the test to be performed that recipient must select to proceed with the test by clicking on "Go".

10

It is worth noting on Fig. 8 that should the recipient experience any difficulties or discomfort during testing, they have the ability to stop the testing at any time by pressing any button or "Stop". This feature is similar to what would occur in a regular clinical visit, where the recipient could indicate to the clinician any issues and the clinician would then stop the testing.

15

Once the recipient has started the system integrity test, the recipient is informed that the test is being performed and is provided with an indication of the progress of the test, see Fig. 9. Upon completion of the test, the recipient is informed that the test has been completed and is asked to proceed to the next stage in the process, see Fig. 10, by clicking "Next".

20

At the completion of the system integrity test, the recipient application undergoes a basic test that checks whether the system integrity test has found any severe system integrity problems. Should such problems be found, the application would inform the recipient of this and suggest that the recipient contact the clinician for further attention. This would then finish the testing process, as the presence of system integrity problems would typically reduce the benefit of obtaining any further test results.

25
30

Following the successful completion of the system integrity test, the recipient is then informed of the next test in the process, namely a psychophysics test, see Fig. 11. Upon commencement of this test, the recipient is given instructions on how the test will be conducted and direction as to how input is required from the recipient. Typically, the input will be in the form of using a mouse to activate the appropriate icon on the computer screen.

35

In the test as shown in Fig. 12, the recipient is asked to press YES when they hear a tone and NO when they do not hear anything. Once this test is complete, the recipient is informed of this and asked to proceed to the next step in the process by pressing "Next" (see Fig. 13).

5

As shown in Fig. 14, when the recipient has completed all tests configured for them, they are informed of this and the application is closed. The recipient presses "End" to close the programme on their computer.

10 As discussed in relation to Fig. 2, upon completion of the configured tests, the results are sent to the clinician for evaluation. This is done by the clinician application which accesses the results file and tests the results in relation to previous or desired results.

15 An example of how this is done is discussed in the process below. This evaluation is related to the results of a T-NRT test performed on each of the recipient's electrodes as set by the clinician.

1. The T-NRT measurements, $T-NRT_n$ for all 22 electrodes of a
20 recipient are stored in the results file following testing by the recipient. These results can be considered as $T-NRT_n(E1)$, $T-NRT_n(E2)$ $T-NRT_n(E22)$.

2. The results of previously measured T-NRT measurements $T-NRT_{n-1}$ are also stored in the clinician application.

25

3. The differences between the current measurements, $T-NRT_n$ and the previously stored measurements, $T-NRT_{n-1}$ are then obtained and compared for each electrode. The clinician would then set a limit as to the difference between the two which is considered acceptable, for example a
30 clinician may consider that a difference of less than 10 current levels is acceptable for all electrodes.

Evaluation:

$|T-NRT_n(E1) - T-NRT_{n-1}(E1)| < 10$ & $|T-NRT_n(E2) - T-NRT_{n-1}(E2)| < 10$
35 $|T-NRT_n(E22) - T-NRT_{n-1}(E22)| < 10$

If the evaluation is positive, the application provides feedback to the clinician that the test is "OK" otherwise the feedback is "NOT OK".

5 As shown in Fig. 15, upon an evaluation of the results as discussed above for all the tests undertaken, the clinician is provided with an interface providing an overview of the results. The overview as shown in Fig. 15 represents the results of the test carried out in Figs. 7-14. As is shown, the test revealed that the recipient's system passed the system integrity test but failed the test of the recipient's threshold levels.

10

The clinician can then decide to view a detailed report of the test, which is shown in Fig. 16. This detailed report shows the measured values obtained in the test against the mapped values, or previously measured values in both a data form and a graphical form. As is evident from this detailed report, the
15 majority of the measured results fall within an acceptable range of the previous results, except at electrode 10, where the difference is quite substantial.

Based upon this result the clinician can decide to proceed in a number of ways. Obviously, if the process indicated that all tests were passed, no specific
20 action would be required and the clinician would inform the recipient of this and ask the recipient to repeat the test in, for example, 6 months time. Should there be some inconsistencies in the results, such as that shown in Fig. 16, the clinician may wish to see the recipient personally or send the recipient further tests to undertake to determine whether the inconsistency is merely an error in
25 recordings or is a concern that requires further monitoring.

As can be appreciated, the present invention provides a system which enables faster, more user friendly and less skill-intensive follow-up care and management of a cochlear implant. This system allows for personal expert
30 care from clinicians and audiologists without the need to onerous clinician visits, and allows the recipient to interact with the clinician from their own home or school.

It will be appreciated by persons skilled in the art that numerous variations and/or modifications may be made to the invention as shown in the specific embodiments without departing from the spirit or scope of the invention
5 as broadly described. The present embodiments are, therefore, to be considered in all respects as illustrative and not restrictive.

Dated this seventh day of November 2002

Cochlear Limited
Patent Attorneys for the Applicant:

F B RICE & CO

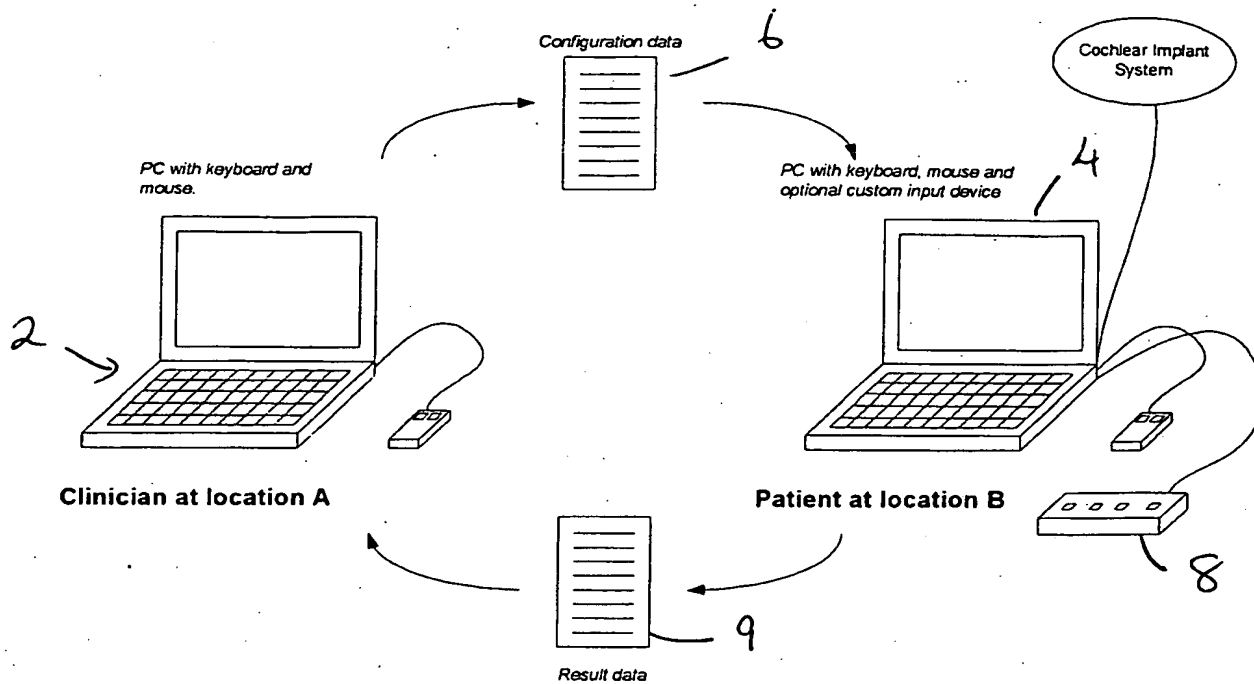


Figure 1

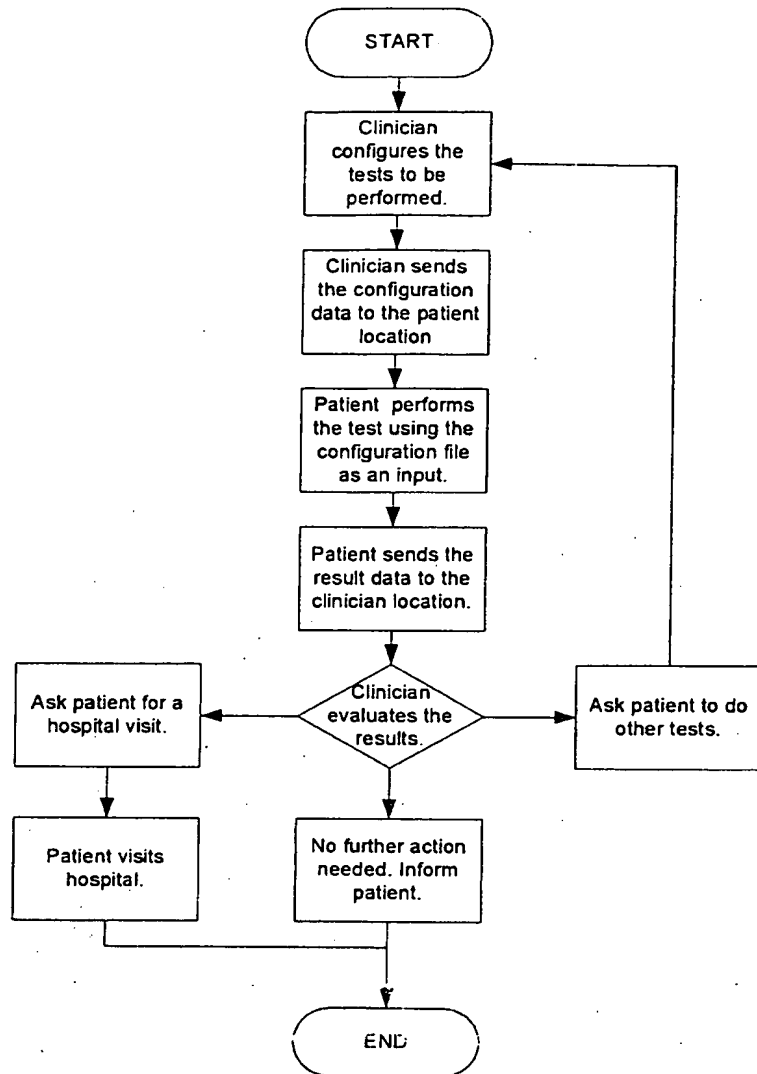


Figure 2

Clinician application

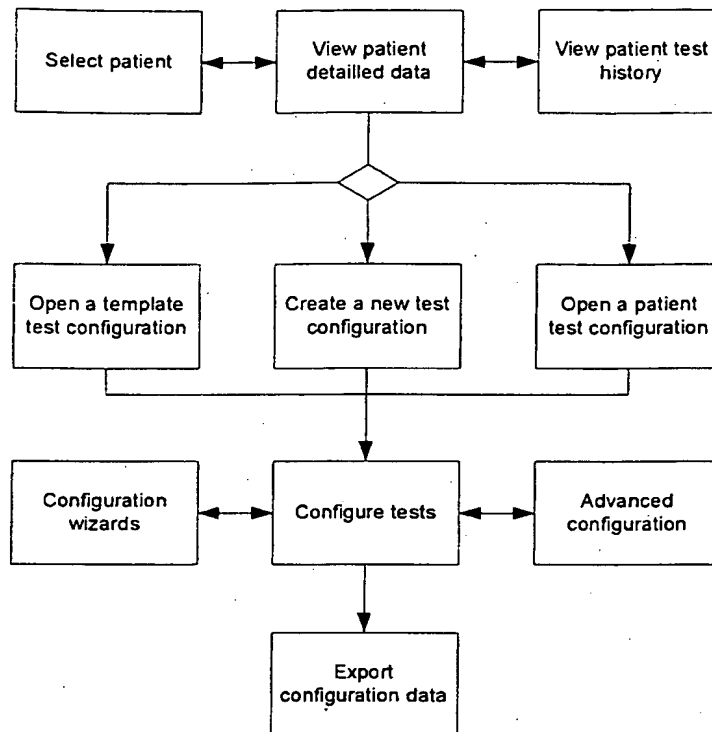


Figure 3

Automated Clinical Assistant Prototype v1.0

1) Recipient:

2) Language:

3) SP Connection:

4) Tests to be performed:

☒ System Integrity

☐ T-NRT

☒ T-levels

☐ C-levels

☐ Loudness balancing

☐ Pitch ranking

☐ Phoneme discrimination

5) Maps to be used:

Save this configuration for remote execution:

Start the recipient's interface:

Automated Clinical Assistant Prototype v1.0

Test results for Bob Peeters

☒ System Integrity: OK

☐ T-NRT

☒ T-levels: NOT OK

☐ C-levels

☐ Loudness balancing

☐ Pitch ranking

☒ Phoneme discrimination

Detailed report:

Redo tests:

Close the application:

Figure 4

Fig. 15

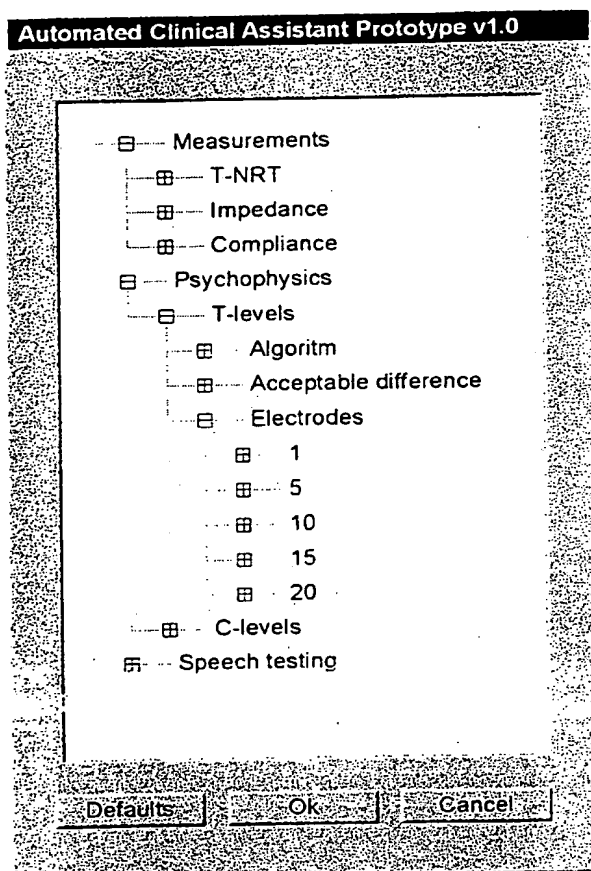


Fig. 5

Patient application

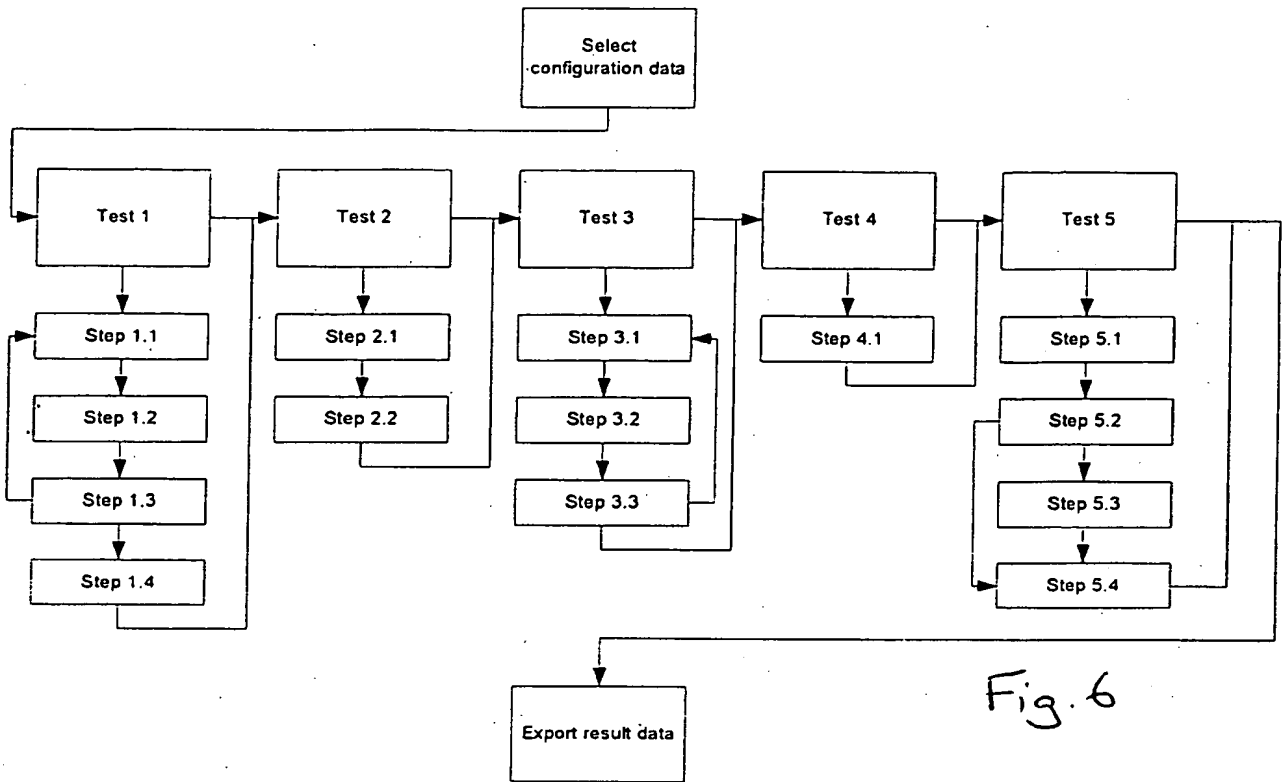


Fig. 6

Welcome, Bob Peeters.

We will lead you through different tests

Ask the audiologist for help if something
is unclear

Press the **NEXT** button on the right
bottom corner to continue

NEXT

Fig. 7

System integrity test

The system will now check your system integrity

Press GO to start the test

Press any button or STOP when you hear something uncomfortable


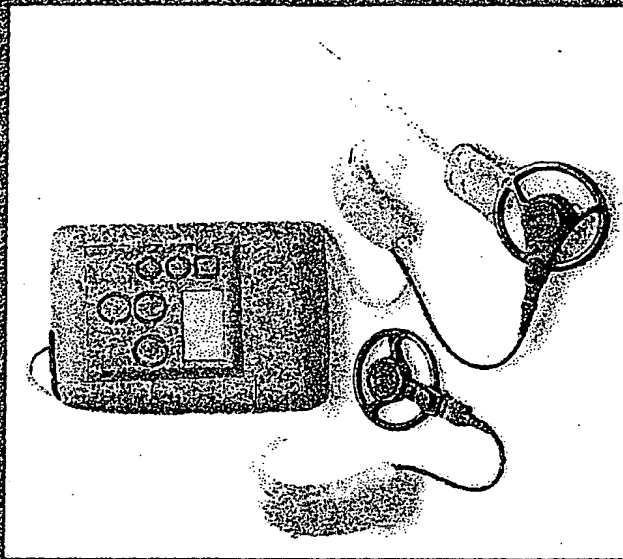


Fig. 8

GO

System integrity test in progress

Press any button or STOP when you hear something uncomfortable

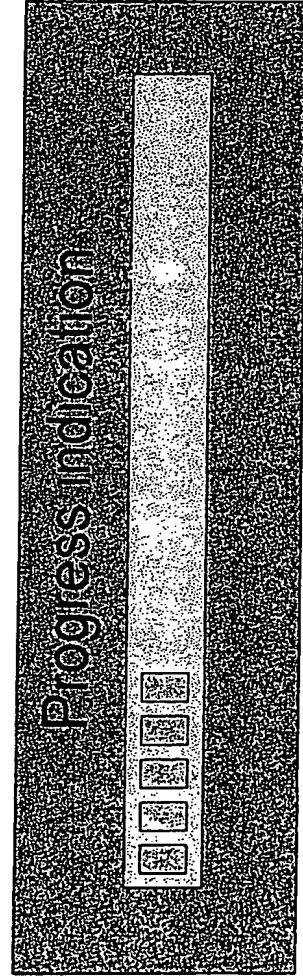

STOP

Fig. 9

System integrity test completed!

Press **NEXT** to continue


NEXT

Fig. 10

Psychophysics test.

The system will now stimulate soft tones.

Press NEXT to continue

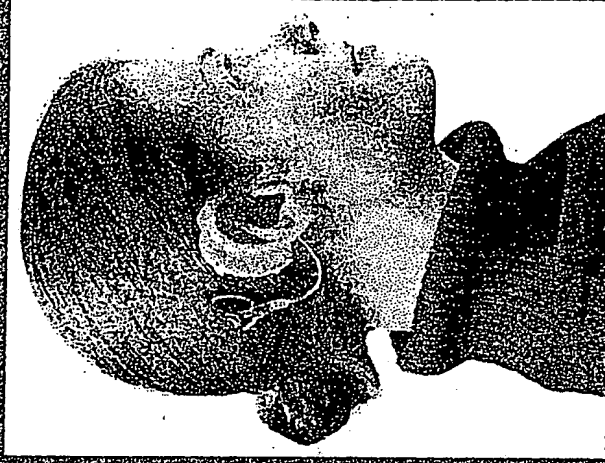

NEXT

Fig. 11

Psychophysics test in progress

Press **PLAY** to start the tone stimulation

Press **YES** when you hear something

Press **NO** when you hear nothing

Press any button or **STOP** when you hear something uncomfortable

You will be unable to hear some tones since we also test at low levels

PLAY

YES

NO

STOP

Fig: 12

Psychophysics test completed!

Press NEXT to continue

NEXT


Fig. 13

All tests completed!

Congratulations, all tests are completed!

Tell your audiologist you have finished

Press END to close the application

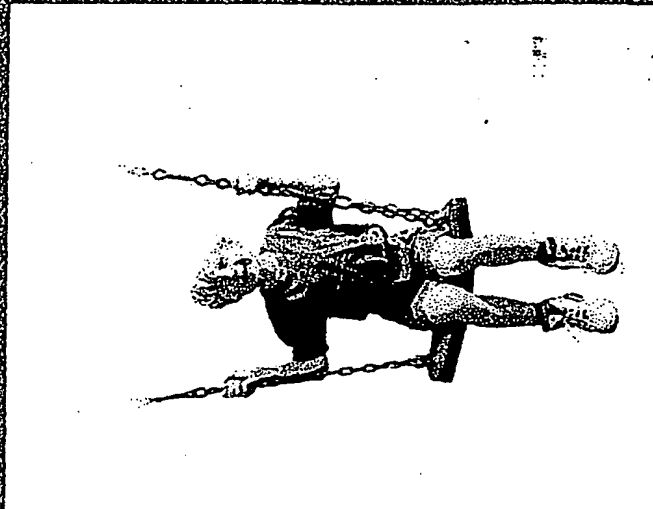

END

Fig. 14

Generated report

Clinician gets a detailed report:

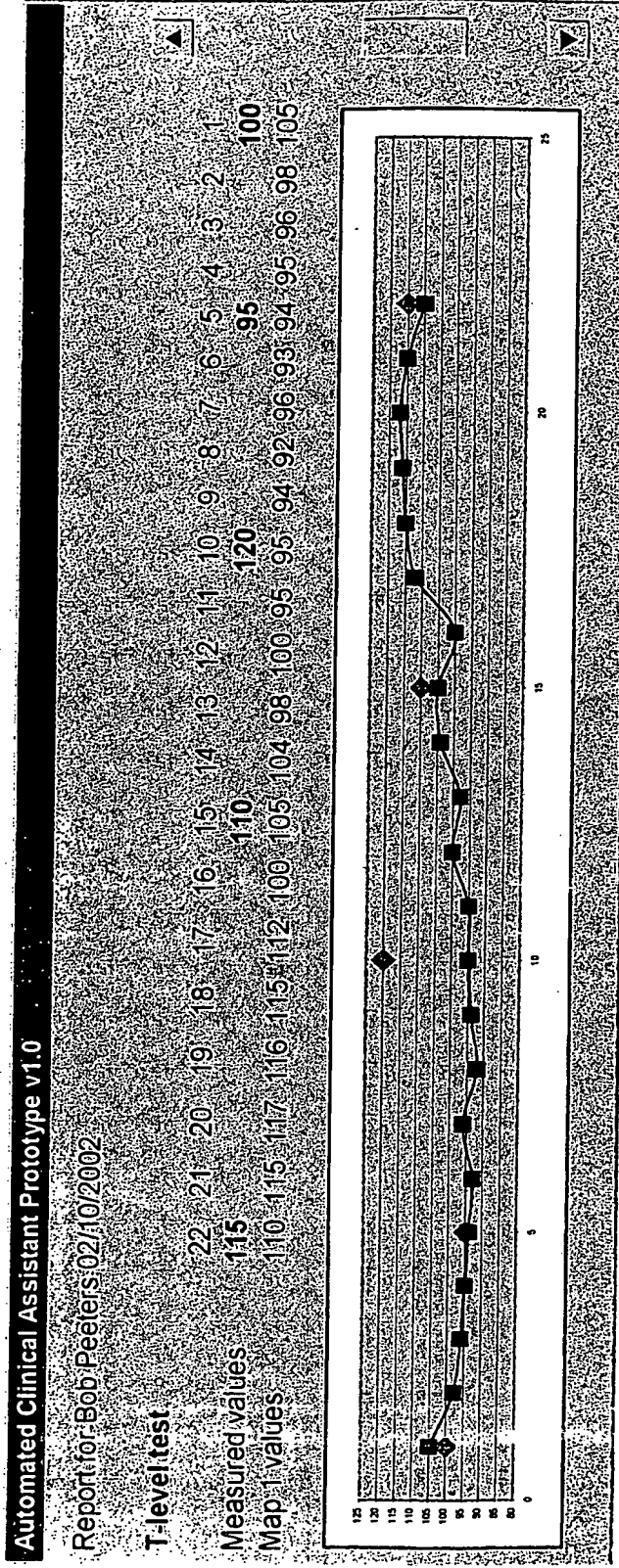


Fig. 16